

# TroClose™ 1200

Novel Integrated  
Port Closure System



Simple  
Secure  
Safe





Cannula

Ribbed cannula



Suture




Absorbable  
anchor

# SIMPLE

# SECURE SAFE

TroClose™ 1200

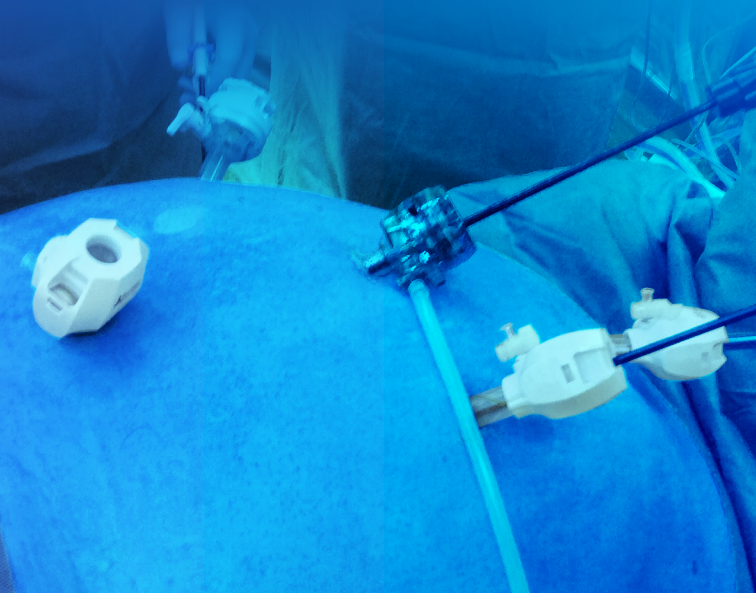

**Novel Integrated  
Port Closure System**

# SECURE SAFE



Obturator

Bladeless  
Obturator tip





# TroClose™ 1200

Novel Integrated Port Closure System

Gordian Surgical's TroClose™1200 is a **novel integrated port closure system** that offers surgeons a simple, secure and safe solution to open and suture closed the abdominal wall during laparoscopic procedures.

Instead of insertion of closure sutures into the fascia at the end of a procedure, Gordian's uniquely designed "two in one" trocar efficiently inserts sutures into the tissue surrounding the cannula at the beginning of the procedure. When the cannula is removed, the surgeon closes the fascia by simply tying the TroClose™1200 sutures together.

[www.GordianSurgical.com](http://www.GordianSurgical.com)

# Technical Specifications

- 12 mm bladeless trocar
- 115 mm ribbed cannula
- Single-use only
- EtO sterilized
- Fully absorbable PLGA anchors
- Braided PGLA absorbable sutures (similar to USP 0)

## Available as

Description	Cat #
TroClose™1200 - Obturator & cannula	GSTC1200
TroClose™1200 - Cannula Only	GSTC1200C



MKT1005, 2020-02-3

**Rx** Prescription Only: U.S. Federal law restricts this device to sale by or on the order of a physician or properly licensed practitioner.

**CE** 2409

**Gordian Surgical Ltd.**  
Misgav Business Park  
M.P. Misgav 2017400, Israel

**Email** [info@gordiansurgical.com](mailto:info@gordiansurgical.com)  
**US** +1-813-240-9448  
**ROW** +972-72-260-7095

**[www.GordianSurgical.com](http://www.GordianSurgical.com)**